


# NEW!

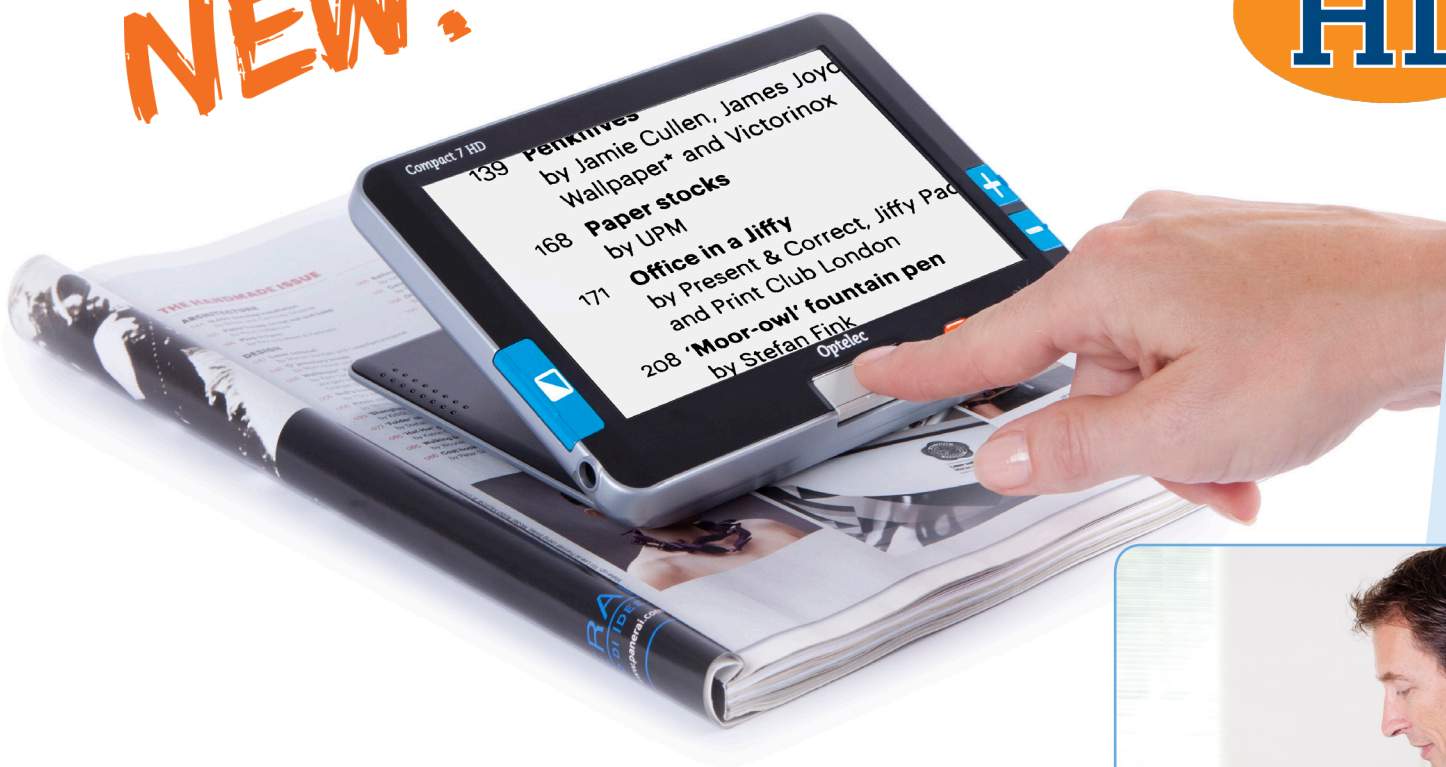

## Optelec Compact 7 HD

The Optelec Compact 7 HD is the newest, truly portable 7-inch widescreen electronic video magnifier with High Definition image quality, extreme simplicity and an attractive style. Integrating the largest screen available in the Optelec Compact product family, it is still small enough to easily carry from one location to another.

*Reading newspapers, letters or viewing photos anywhere at home, in the office, at school or while on vacation is easier than ever!*

**Optelec**  
Life is worth enjoying

# Large screen in a compact design



The continuous magnification range starting from 2X allows you to see the smallest details of photos or objects on the ergonomically positioned, high contrast 7-inch widescreen. Maximize your reading ability with its increased field of view.

The Optelec Compact 7 HD incorporates dual viewing modes for reading books and long documents when raised open, or closed for distance needs such as organizing around the house or viewing a wall temperature display.



## Read with ease & quality immediately

With only "one-touch" to open, the Compact 7HD tilts to a raised position for a comfortable and ergonomic viewing angle. To use, simply place it on the object you want to view and press the large button underneath the screen. It will automatically magnify text or images with your preferred settings. For an effortless user experience, all Compact 7 HD controls are easy to see and operate.

The 8 megapixel High Definition and Auto-Focus camera delivers superior image quality and keeps the text or photos on the screen sharp and clear. A patented LED lighting system eliminates glare and light reflections when reading from glossy documents and displays.





# Take it anywhere

Unlike desktop video magnifiers, the Optelec Compact 7 HD allows you to comfortably read letters and books on your couch or in your favorite chair. Weighing less than 1.4 pounds, the size and travel case of the Compact 7 HD make it easy to take anywhere.



## Customizable usability

Adjust the settings to meet your specific needs in the special large icon menu. Use the intuitive, easy to operate buttons to select the brightness, power options or check the remaining battery percentage. Save up to five preferred viewing modes including full color photo and four of the 18 different high contrast color combinations available.






High Definition Camera



Auto-Focus

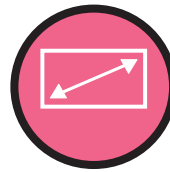

7-inch widescreen



Reading



Ergonomic  
screen position




2X-24X  
magnification



High contrast  
viewing modes



Large icon  
menu

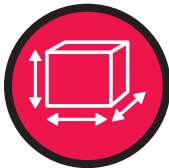

7.5 x 5.1 x 1.1 in  
190 x 130 x 26.7 mm



Average  
4 hours use



22.6 ounces  
1.41 pounds



Snapshot



*Connect with us!*



*Available through our Authorized Optelec Dealer Network:*

**Optelec U.S. Inc.**

**Tel: 800.826.4200**

**Fax: 800.368.4200 E: [info@optelec.com](mailto:info@optelec.com)**

**[www.Optelec.com](http://www.Optelec.com)**

**Optelec U.S. Inc.**

3030 Enterprise Court, Suite C - Vista, CA 92081

**Optelec has offices in:**

THE NETHERLANDS – GERMANY – BELGIUM – UNITED KINGDOM – USA – CANADA

*“Improving the quality of life for the visually impaired  
– reaching out with simple and effective solutions.”*